# Product Characterization Sheet HE3033





#### Classification

Plateability	Suspension	
Number of days plateable	NA	
Confluency	NA	
P450 Inducibility	NA	
Number of donors	1	
Enterocyte Size	15.2 ± 3.1 microns	

## Post-thaw Viability and Yield

Viability	85 %	
Yield	3.7 million	

#### **Donor Demographics**

Candan	Mala	
Gender	Male	
Age	32 years	
Race	Hispanic	
Cause of death	Anoxia	
BMI	21.0	
Smoking	No	
Alcohol	No	
Substance abuse	Yes	
Medical history	Diabetes, HTN	
Infectious diseases	HBV-, HCV-, HIV-, CMV+	


<u>Characterization:</u> Enterocytes were thawed using 37°C CERM<sup>TM</sup> and centrifuged for 10 minutes at 100g. After removing the supernatant, enterocytes were re-suspended in 4°C HQM<sup>TM</sup> and counted for viability and yield using the Trypan Blue exclusion method.

### **Drug Metabolism Activity**

Drug Metabolizing Enzyme	Substrate (μΜ)	Incubation Time (minutes)	Metabolite Quantified	Activity (pmol/minute/million cells)
CYP3A4	Midazolam (20)	120	1-Hydroxymidazolam	1.03 ± 0.09
	Testosterone (200)	120	6β-Hydroxytestosterone	28.4 ± 7
CYP2J2	Astemizole (50)	120	O-Desmethyl astemizole	0.7 ± 0.05
CES2	Irinotecan (50)	120	SN38	$0.11 \pm 0.03$
CYP2C9	Diclofenac (25)	120	4-Hydroxydiclofenac	2.11 ± 0.05
CYP2C19	S-Mephenytoin (250)	120	4'-Hydroxymephenytoin	1.13 ± 0.13
Sulfotransferase	7-Hydroxycoumarin (100)	120	7-Hydroxycoumarin sulfate	14.2 ± 1.1
UGT	7-Hydroxycoumarin (100)	120	7-Hydroxycoumarin glucuronide	29.6 ± 2.5

<u>CYP450 Activity Assessment:</u> The enterocytes were incubated at a cell density of 1.5 million enterocytes/mL in a 96-well plate (150,000 enterocytes/well) for the designated time duration of 120 minutes with isoform-selective substrates. The metabolites were identified and analyzed using LC-MS/MS.

## Photomicrographs (100X, Brightfield)

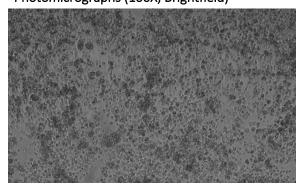

<u>Enterocyte Morphology:</u> Enterocyte cell size is approximately 12 – 17 microns in diameter. Lymphocytes are approximately 4 microns and may also be observed in the enterocyte cell population.

IVAL cell culture media and tissue culture plates used in this evaluation:

- Recovery of thawed enterocytes Cat. No. 81056 CERM™ Cryopreserved Enterocyte Recovery Media, 50 mL tube
- Suspension and incubation of enterocytes Cat. No. 81039 HQM™ Hepatocyte and Enterocyte Incubation Media, 50 mL tube

To inquire about our products and services or for technical questions please contact:

• In Vitro ADMET Laboratories by phone at +1 (866) 458-1094 or +1 (410) 869-9037 or email at info@invitroadmet.com