# Product Characterization Sheet HE3019





### Classification

Plateability	Suspension	
Number of days plateable	NA	
Confluency	NA	
P450 Inducibility	NA	
Number of donors	1	
Enterocyte Size	15.4 ± 2.0 microns	

## Post-thaw Viability and Yield

Viability	89 %	
Yield	0.7 million	

#### **Donor Demographics**

Male		
61 years		
Caucasian		
Anoxia		
30.0		
No		
No		
No		
Melanoma		
HBV-, HCV-, HIV-, CMV+, EBV (IgG)+		


<u>Characterization:</u> Enterocytes were thawed using 37°C CERM<sup>TM</sup> and centrifuged for 10 minutes at 100g. After removing the supernatant, enterocytes were re-suspended in 4°C HQM<sup>TM</sup> and counted for viability and yield using the Trypan Blue exclusion method.

## **Drug Metabolism Activity**

Drug Metabolizing Enzyme	Substrate (μΜ)	Incubation Time (minutes)	Metabolite Quantified	Activity (pmol/minute/million cells)
CYP3A4	Midazolam (20)	120	1-Hydroxymidazolam	0.30
	Testosterone (200)	120	6β-Hydroxytestosterone	5.27 ± 1.34
CYP2J2	Astemizole (50)	120	O-Desmethyl astemizole	0.26
CES2	Irinotecan (50)	120	SN38	0.30
CYP2C9	Diclofenac (25)	120	4-Hydroxydiclofenac	0.24
CYP2C19	S-Mephenytoin (250)	120	4'-Hydroxymephenytoin	0.11
Sulfotransferase	7-Hydroxycoumarin (100)	120	7-Hydroxycoumarin sulfate	1.89
UGT	7-Hydroxycoumarin (100)	120	7-Hydroxycoumarin glucuronide	1.47

<u>CYP450 Activity Assessment:</u> The enterocytes were incubated at a cell density of 1.5 million enterocytes/mL in a 96-well plate (150,000 enterocytes/well) for the designated time duration of 120 minutes with isoform-selective substrates. The metabolites were identified and analyzed using LC-MS/MS.

# Photomicrographs (100X, Brightfield)

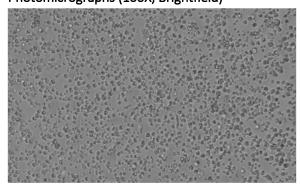

<u>Enterocyte Morphology:</u> Enterocyte cell size is approximately 12 – 17 microns in diameter. Lymphocytes are approximately 4 microns and may also be observed in the enterocyte cell population.

IVAL cell culture media and tissue culture plates used in this evaluation:

- Recovery of thawed enterocytes Cat. No. 81056 CERM™ Cryopreserved Enterocyte Recovery Media, 50 mL tube
- Suspension and incubation of enterocytes Cat. No. 81039 HQM™ Hepatocyte and Enterocyte Incubation Media, 50 mL tube

To inquire about our products and services or for technical questions please contact:

• In Vitro ADMET Laboratories by phone at +1 (866) 458-1094 or +1 (410) 869-9037 or email at info@invitroadmet.com